# Tolaccea

## PRODUCT MANUAL

www.tolaccea.com

### TOLACCEA ANKLE ICE PACK PRO

Hot and  
Cold Therapy

Natural and  
Expert-grade Gel

Soft Lycra  
Fabric

Fit for Left &  
Right Ankles

#### WARNING:


To avoid danger of suffocation, keep this away from children and babies!

Company Name: REPLAB SARL  
Address and Post Code: 99 AVENUE ACHILLE  
PERETTI 92200 NEUILLY-SUR-SEINE France  
Contact: XIAOLAN LIU  
Phone: +33767189223  
Email : EUREP@REPLAB.CC

Company Name: REPLAB CONSULTING LTD  
Address and Post Code: 376 NEWPORT ROAD CARDIFF,  
WALES, CF23 9AE  
Contact: HUIHUA LU  
Phone: +447395577022  
Email: UKREP@REPLAB.CC
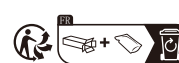
## EN

### FUNCTIONS

1. Tolaccea hot&cold therapy ice sleeve provides 360°coverage. Our ColdFocus Technology provides a double layer of Cold Therapy, ensuring extra cooling precisely for the ankle.
2. The flexible ice foot sleeve is in ergonomic design, perfect for foot and it is wearable without straps.
3. This ice foot sleeve is reusable, also microwave and freezer safe, allowing you to use hot or cold therapy depending on your needs.
4. The ice foot sleeve is made of non-toxic material, latex-free, and leak proof, with no smells after use.
5. It can help to relieve pain for : achilles tendonitis, plantar fasciitis, sprained ankles or feet.

### OPERATION INSTRUCTION

#### Cold therapy

- Place the foot sleeve flat in a refrigerator at -20 °C for 1-2 hours for immediate use. If colder temperature needed, freeze it longer with lower temperature.
- Take out the foot sleeve from freezer, wipe off the moisture and pull it over the required area.
- Wear foot sleeve and keep in place for approximately 15-20

#### Hot therapy

- Start with room temperature pack, place the foot sleeve flat in microwave with a clean plate.
- Microwave the foot sleeve for 10 seconds, if more heat is needed microwave in 5 second increments.
- Wear it and keep in place for approximately 15-20 minutes.

### HOW TO CLEAN IT

Hand wash with damp towel. Do not submerge sleeve in water.  
Lay flat and air dry, do not machine wash.

### WARNING!

1. Limit cold/hot therapy to 30 minutes to prevent frostbite or burns. If not cold enough, refreeze at -20°C.
2. Do not microwave over 45 seconds. Check temperature before use; cool if too hot to avoid burns.
3. Repeated heating may dry out the gel, reducing lifespan. Monitor while microwaving to avoid overheating. Air dry if damp—never use a dryer.
4. For external use only. Avoid areas with poor circulation, open wounds, or infections unless advised by a professional.
5. If gel contacts skin, wash thoroughly. Seek medical advice for injuries, and discontinue if discomfort persists.
6. Store in therapy wrap when not in use to prevent condensation. Ensure dry before storage.
7. Keep away from children and pets.

### AFTER SALE

Tolaccea guarantees 1 year warranty for the quality problem. Any questions please freely email us: service@tolaccea.com.

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